# STORZ

## KARL STORZ — ENDOSKOPE

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**1 Important information for  
users of KARL STORZ  
devices**

**It is recommended to check the suitability of  
the product for the intended procedure prior  
to use.**

Thank you for placing your trust in the name  
KARL STORZ. Like all of our other products, this  
product is the result of years of experience and  
great care in manufacture. You and your organiza-  
tion have decided in favor of a modern, high quali-  
ty piece of equipment from KARL STORZ.

This instruction manual is intended to aid you  
in the proper installation, connection, and  
operation of this device. All required details and  
all actions on your part are clearly presented  
and explained. Please read this manual carefully  
before proceeding with the equipment. Keep this  
instruction manual available for ready reference in  
a convenient location near the equipment.

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

The IMAGE1 S comprises IMAGE1 S CONNECT™ (TC 200) and a minimum of one and maximum of three additional LINK modules. The IMAGE1 S H3-LINK (TC 300), IMAGE1 S X-LINK (TC 301) and IMAGE1 S D3-LINK™ module (TC 302) are available as LINK modules.

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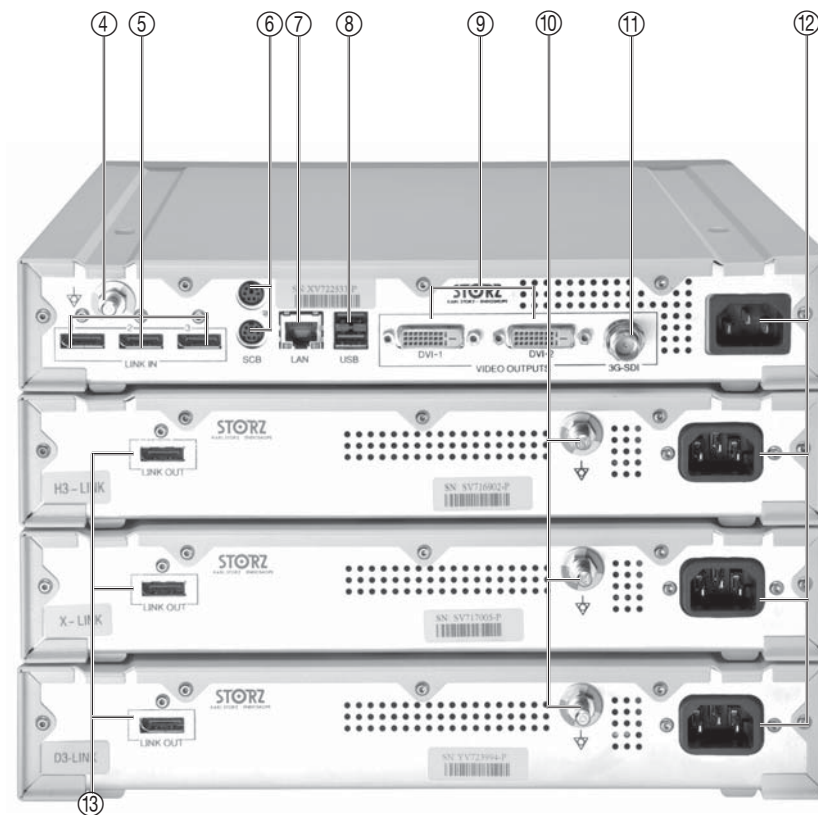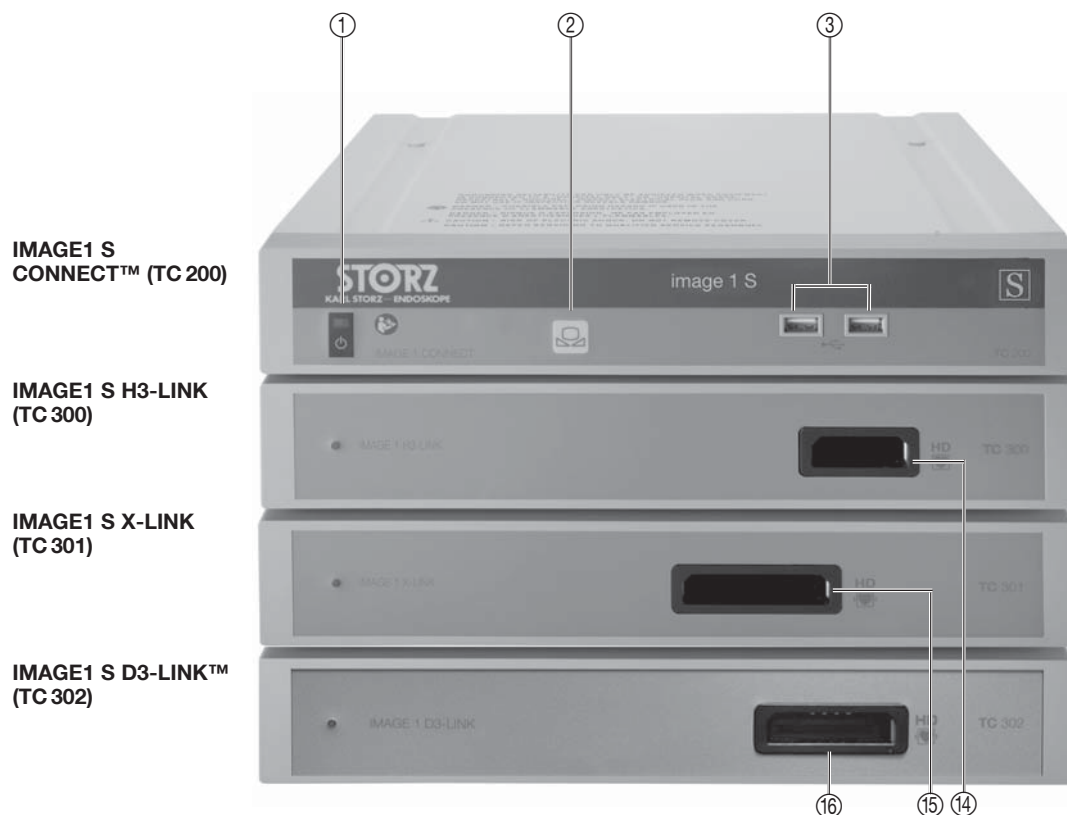**3 Controls, displays,  
connectors,  
and their uses****IMAGE1 S CONNECT™ in combination with  
the H3-LINK, X-LINK and D3-LINK™ modules**

- ① Power switch
- ② White balance
- ③ USB port
- ④ Potential equalization connector
- ⑤ LINK inputs
- ⑥ SCB input/output sockets  
(KARL STORZ Communication Bus)
- ⑦ LAN interface
- ⑧ USB interface
- ⑨ HD DVI-D video output sockets
- ⑩ Potential equalization connector
- ⑪ 3G-SDI video output
- ⑫ Power cord socket
- ⑬ LINK outputs
- ⑭ H3-LINK camera connector
- ⑮ X-LINK camera connector
- ⑯ D3-LINK™ camera connector

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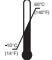
<b>1</b>	<b>Important information for users of KARL STORZ devices</b> .....	<b>III</b>
<b>2</b>	<b>Images of the equipment</b> .....	<b>IV</b>
2.1	Possible combinations for IMAGE1 S .....	IV
<b>3</b>	<b>Controls, displays, connectors, and their uses</b> .....	<b>V</b>
<b>4</b>	<b>Symbols employed</b> .....	<b>VI</b>
4.1	Symbols for operation.....	VI
4.2	Label and packaging symbols.....	VII
<b>5</b>	<b>General information</b> .....	<b>4</b>
5.1	Description of the devices.....	4
5.2	Property rights .....	4
<b>6</b>	<b>Safety instructions</b> .....	<b>5</b>
6.1	Explanation of warnings and cautions.....	5
6.2	General safety information .....	6
6.3	Indications for use .....	10
6.3.1	Intended use.....	10
6.4	Contraindications.....	10
6.5	Patient profile.....	11
6.6	User profile (specialist physician) .....	11
6.7	User profile (assistant) .....	12
6.8	Safety precautions at the site of installation .....	13
6.9	Safety precautions when operating the device.....	13
<b>7</b>	<b>Initial operation</b> .....	<b>14</b>
7.1	Basic equipment.....	14
7.1.1	TC 200 IMAGE1 S CONNECT™.....	14
7.1.2	TC 300 IMAGE1 S H3-LINK.....	14
7.1.3	TC 301 IMAGE1 S X-LINK .....	14
7.1.4	TC 302 IMAGE1 S D3-LINK™ .....	14
7.2	Unpacking the equipment.....	14
7.3	Preparing the device for use .....	15
7.3.1	Installation.....	15
7.3.2	Connecting the potential equalization.....	15
7.3.3	Connecting the IMAGE1 S LINK module with IMAGE1 S CONNECT™ .....	15
7.3.4	Connecting SCB.....	15

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7.3.5	Connecting the monitor .....	16
7.3.6	Connecting USB devices.....	16
7.3.7	Connecting the power supply .....	16
7.3.8	Switching on the IMAGE1 S CONNECT™ ..	17
7.3.9	Connecting USB storage media .....	17
7.3.10	Connecting the camera head/ video endoscope .....	17
7.3.11	Performing the white balance .....	17
7.3.12	Test for proper functioning .....	18
7.3.13	Installation overview .....	19
<b>8</b>	<b>Operating instructions.....</b>	<b>20</b>
8.1	Basic information on operation .....	20
8.1.1	Keyboard and mouse .....	20
8.1.2	Camera head buttons.....	20
8.2	Dashboard.....	20
8.2.1	Description of functions of the symbols used .....	21
8.3	Live Menu.....	25
8.3.1	Description of functions of the symbols used .....	26
8.4	Setup Menu.....	31
8.4.1	General Settings (IMAGE1 S CONNECT™) .....	31
8.4.2	Presets .....	35
8.4.3	LINK-specific settings .....	35
<b>9</b>	<b>Cleaning, disinfection, care and sterilization .....</b>	<b>36</b>
9.1	General safety information .....	36
9.2	Wipe-down disinfection of devices.....	37
<b>10</b>	<b>Service and repair.....</b>	<b>38</b>
10.1	Maintenance and safety check .....	38
10.1.1	Maintenance .....	38
10.1.2	Safety check.....	38
10.2	Servicing and repair .....	39
10.3	Disposal.....	39
10.4	Repair program.....	40
10.5	Limitation of liability.....	41
10.6	Warranty .....	41
<b>11</b>	<b>Technical description.....</b>	<b>42</b>
11.1	Technical data.....	42
11.2	Technical documentation .....	43
11.3	Standard compliance.....	43

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11.4	Directive compliance for IMAGE1 S (TC 200, TC 300, TC 301 and TC 302) .....	44
11.5	Software ownership and licensing .....	44
<b>12</b>	<b>Spare parts, recommended accessories.....</b>	<b>46</b>
12.1	Compatible camera heads for IMAGE1 S H3-LINK TC 300 with S technologies .....	46
12.2	Compatible camera heads for IMAGE1 S H3-LINK TC 300 without S technologies .....	46
12.3	Compatible camera heads for IMAGE1 S X-LINK TC 301 with S technologies.....	47
12.4	Compatible camera heads for IMAGE1 S X-LINK TC 301 without S technologies.....	47
12.5	Compatible video endoscopes for IMAGE1 S X-LINK TC 301 with S technologies .....	47
12.6	Compatible video endoscopes and camera heads for IMAGE1 S D3-LINK™ TC 302 with S technology .....	49
12.7	Recommended accessories .....	50
12.8	Operable SCB light sources.....	51
12.9	Operable SCB insufflators.....	51
12.10	Compatible 3D monitors.....	51
<b>13</b>	<b>Information on electromagnetic compatibility (EMC).....</b>	<b>52</b>
13.1	EMV tables for IMAGE1 S.....	54
<b>14</b>	<b>Subsidiaries .....</b>	<b>66</b>

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## 5 General information

### 5.1 Description of the devices

The IMAGE1 S takes a unique and innovative approach to system conception. The modular design principle allows the user to create a system which is tailored to his needs. The IMAGE1 S CONNECT™ module is at the heart of the IMAGE1 S and acts as a communication interface between the image and documentation devices.

For use, one IMAGE1 S CONNECT™ plus at least one LINK module is required which receives and forwards the signal. The Light Source Control and S technologies for image optimization provide users with even more comprehensive support. An intuitive user interface assures simple and quick access as well as rapid understanding of menu navigation.

The system offers numerous state-of-the-art functions. These include:

- IMAGE1 S visualization modes to assist the user via contrast enhancement and/or more even brightness perception
- Light Source Control for automatic adjustment of the light source
- Integrated image and video capture for simple documentation

Functions are described in greater detail in § 8.4.

### 5.2 Property rights

Protected by at least one of the following US Patents: 7,821,530; 8,274,559; 7,520,853; 5,716,323; 5,788,688; 5,913,817; 6,397,286; 6,484,221; 7,212,227; US 8,089,509; US 8,199,188. For more information, visit: [www.karlstorz.com](http://www.karlstorz.com).

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**6 Safety instructions****6.1 Explanation of warnings  
and cautions**

Safety instructions are measures intended to protect the user and patients from the risks which could arise through the use of the system.

Please read this manual and follow the instructions carefully. The words **Warning**, **Caution** and **Note** convey special meanings. Wherever they are used in this manual, the accompanying text should be carefully reviewed to ensure the safe and effective operation of the system. To make these words stand out more clearly, they are accompanied by a pictogram.

**WARNING:** A Warning indicates that the personal safety of the patient or physician is at risk.

**CAUTION:** A Caution indicates that particular service procedures or precautions must be followed to avoid possible damage to the device.

**NOTE:** A Note indicates special information about operating the device or clarifies important information.

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**6.2 General safety information**

**WARNING:** Risk of injury and damage to the products: Failure to observe and follow this instruction manual and all of the instruction manuals of the products used in combination can result in injury to patients, users and third parties as well as damage to the product. Please read all relevant instruction manuals carefully and always follow the instructions given precisely. Check the functioning of the products used in combination.

**WARNING:** Test this device prior to each surgical procedure to ensure that it functions correctly. In the event that the image becomes unusable during surgery, the camera may be disengaged from the endoscope and the procedure continued optically. If this is not possible, it is up to the discretion of the surgeon how best to proceed. Availability of a spare system is recommended.

**WARNING:** Danger of an electric shock! This device may only be connected to a socket with protective earthing. To reduce the risk of electric shock, keep the device out of reach of patients.

**WARNING:** Grounding reliability can only be achieved when the equipment is connected to a 'Hospital Only' or 'Hospital Grade' outlet (i.e., approved for use in an operating room environment). The plug and cord should be routinely inspected. Do not use if damaged.

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 **WARNING:** Consult and follow the storage, transport and operating conditions prior to storing, transporting or operating, respectively. If the device is stored in a warm, humid environment, ensure that it has suitably equilibrated to the operating environment prior to use.

 **WARNING:** To reduce the risk of electrical shock, do not remove cover of unit. Refer servicing to qualified personnel. Removal of covers by unauthorized personnel will void the equipment's warranty.

 **WARNING:** Perceived electrical shocks may be experienced when using the camera system in conjunction with electrosurgical units.

 **WARNING:** High energy radiated light through endoscopes may give rise to high temperatures in front of the light outlet and to the tip of the endoscope. To minimize the risk of burns, always adjust the light source to the minimum illumination intensity necessary to achieve optimum illumination of the endoscopic scene when coupled to the video camera.

 **WARNING:** Refer to the appropriate section of this manual for validated cleaning and disinfection instructions.

 **CAUTION:** Do not store liquids on or above the unit.

 **CAUTION:** In the USA, federal law restricts this device to sale or use by or on the order of a physician.

 **NOTE:** Disposal of these system products, at the end of their useful life, must be in accordance with local regulations.

 **NOTE:** Do not discard as unsorted municipal waste.

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 **NOTE:** Discard as electrical/electronic waste; recycle or reuse accordingly.

 **WARNING:** No modification of this equipment is allowed.

 **NOTE:** Consult local authorities for reuse/recycle instructions.

 **NOTE:** The presence of maximum mains voltage on SIP/SOPs is not considered a normal condition.

 **WARNING:** Accessory equipment connected to the analog and digital interfaces must be certified according to the respective IEC standards (e.g. IEC 60950 for data processing equipment and Clause 16 of 60601-1 3rd edition for medical equipment). Furthermore, all configurations shall comply with the system standard Clause 16 of 60601-1 3rd edition. Any person who connects additional equipment to the signal input part or signal output part configures a medical system, and is therefore responsible for ensuring that the system complies with the requirements of the system standard Clause 16 of 60601-1 3rd edition. If in doubt, consult the technical service department or your local representative.

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**WARNING:** To ensure safe operation, do not simultaneously touch the device output connectors and the patient.

**WARNING:** When endoscopes are used with energized endoscopically used accessories, the patient leakage currents may be additive. This is particularly relevant when using type CF devices. In this case, a Type CF endoscopically used accessory should be used in order to minimize the total patient leakage current.

**WARNING:** During discharge of a defibrillator, the image displayed may take up to 3 seconds to recover and in some cases a power cycle may be required to restore the image.

**WARNING:** The device should be placed on a flat surface or in a video cart. Do not block air inlets/outlets as overheating may result in fire hazard or device failure during procedure.

**WARNING:** Failure to perform proper white balance of system may result in color distortion.

**WARNING:** The operator is responsible to limit physical access to the IMAGE1 S from unauthorized persons.

**WARNING:** In the case of loss of surgical image for greater than 3 seconds, a power cycle may be required to restore the image.

**WARNING:** When using high frequency surgical equipment, keep the working part of the active electrode in the field of view of the operator to avoid accidental high frequency burns. Avoid contact between the active electrode and metal parts of the endoscope.

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**6.3 Indications for use****6.3.1 Intended use**

IMAGE1 S is a camera control unit (CCU) for use with camera heads or video endoscopes for the visualization and documentation of endoscopic and microscopic procedures.

**6.4 Contraindications**

No contraindications relating directly to the medical device are currently known. The responsible physician must decide whether the foreseen application is admissible based on the general condition of the patient.

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**6.5 Patient profile**

Use of the IMAGE1 S is not limited to a certain patient profile (sex, age, weight etc.). The attending physician must assess a patient's health with regard to the endoscopic intervention and declare it adequate.

**6.6 User profile (specialist physician)**

- Recognized medical skills for endoscopic interventions in the relevant specialization
- Adequate powers of comprehension to assess the situation rationally at all times during the application
- Adequate language skills in at least one of the languages used on the device and in the instruction manual
- Be thoroughly trained in the operation and use of the device
- Knowledge of the contents of the instruction manual
- No physical impairments that could diminish perception of activation and visual message

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
**6.7 User profile  
(assistant)**

- Completed training as medical assistant
- Adequate language skills in at least one of the languages used on the device and in the instruction manual
- Be thoroughly trained in the operation and use of the device
- Knowledge of the contents of the instruction manual
- No physical impairments that could diminish perception of activation and visual message

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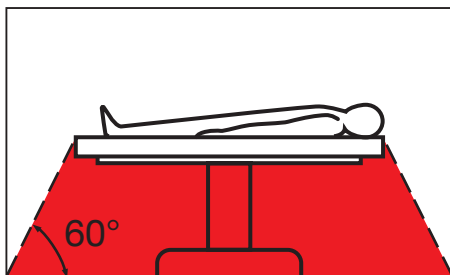### 6.8 Safety precautions at the site of installation

The unit may only be used in medical rooms whose electrical systems have been installed in accordance with applicable national regulations. The system is not intended for use in areas where there is a risk of explosion. This means, for example:

Wherever easily combustible and explosive inhalation anesthetics and their mixtures are used, the device may not be operated in the demarcated hazard zone. This also applies to easily combustible and potentially explosive chemicals, e.g., skin disinfectants and fast-acting surface disinfectants.

The unit is equipped with a connector for potential equalization. It should be connected in accordance with the applicable national regulations.

### 6.9 Safety precautions when operating the device

It is the user's responsibility to make sure the device is safe and operating properly before using the device.

During treatment with the equipment, the patient must be treated and kept under observation with the usual medical care. This includes keeping a check on the progress of treatment, as well as monitoring the vital levels and the anesthetic.

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### 7.1 Basic equipment

#### 7.1.1 TC 200 IMAGE1 S CONNECT™

- 1 power cord, length 300 cm
- 1 DVI connecting cable, length 300 cm
- 1 SCB connecting cable, length 100 cm
- 1 3G-SDI Cable\*
- 1 KARL STORZ USB stick
- 1 Keyboard with touchpad
- 1 instruction manual

#### 7.1.2 TC 300 IMAGE1 S H3-LINK

- 1 power cord, length 300 cm
- 1 LINK cable, 20 cm

#### 7.1.3 TC 301 IMAGE1 S X-LINK

- 1 power cord, length 300 cm
- 1 LINK cable, 20 cm

#### 7.1.4 TC 302 IMAGE1 S D3-LINK™

- 1 power cord, length 300 cm
- 1 LINK cable, 20 cm

## 7.2 Unpacking the equipment

Carefully remove the device and accessories from the packaging. Check for missing items and evidence of damage during shipment. Please file any complaints with the manufacturer or supplier immediately.

If possible, retain the original packing materials for later use; these may be useful if the unit has to be transported.

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### 7.3 Preparing the device for use

#### 7.3.1 Installation

**i NOTE:** Up to three IMAGE1 S LINK modules can be connected to the IMAGE1 S CONNECT™.

**i NOTE:** Do not cover the lateral air vents. Maintain a minimum distance of 6 cm to other devices.

Place the IMAGE1 S CONNECT™ and the corresponding IMAGE1 S LINK modules on a flat surface.

#### 7.3.2 Connecting the potential equalization

**i NOTE:** The unit's ground line should be installed by a qualified electrician.

The IMAGE1 S CONNECT™ and the IMAGE1 S LINK modules are equipped with connectors for potential equalization (see fig. A and B).

1. Connect the potential equalization line to the connectors for potential equalization on all modules.
2. Connect the potential equalization line to the connector in the treatment room.

#### 7.3.3 Connecting the IMAGE1 S LINK module with IMAGE1 S CONNECT™

Using the LINK cable supplied, connect the LINK module to the IMAGE1 S CONNECT™ (see fig. C).

**i NOTE:** The IMAGE1 S CONNECT™ inputs are numbered. Please connect the LINK modules in ascending order.

#### 7.3.4 Connecting SCB

Pull back the protection device of the SCB connector and insert the connector into one of the SCB sockets of the TC 200 IMAGE1 S CONNECT™. Connect the other end of the cable to the KARL STORZ-SCB (KARL STORZ Communication Bus) control unit or other SCB units (see KARL STORZ-SCB System Instruction Manual) (see fig. D).

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**7.3.5 Connecting the monitor**

Connect the monitor with the IMAGE1 S CONNECT™ (see fig. E). Depending on the types of connectors provided on the monitor, this connection can be generated using either:

- DVI-D output
- 3G-SDI output

The primary working screen should be connected to DVI1.

**i NOTE:** The connected monitor must support a resolution of 1920 x 1080 pixels.

**i NOTE:** A compatible 3D monitor must be connected when using a D3-LINK™ in conjunction with 3D video endoscopes (see chapter 12.10).

**7.3.6 Connecting USB devices**

The USB ports of the IMAGE1 S CONNECT™ are suitable for connecting up the following peripherals (see fig. F):

- Keyboard/Mouse
- Printer
- Footswitch
- USB to ACC adaptor
- USB flash drives, FAT32

**7.3.7 Connecting the power supply**

Before connecting the devices to the electrical supply line, check that the values on the power sockets concur with the values on the devices' identification plates.

1. Insert the device plug of the power cord into the power supply cord receptacle of the IMAGE1 S CONNECT™ (see fig. G).
2. Insert the device plug of the power cord into the power supply cord receptacle of the IMAGE1 S LINK modules (see fig. H).
3. Insert the modules' power cords into the power sockets.

**! WARNING:** Always use a hospital grade power cord with this camera.

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**7.3.8 Switching on the IMAGE1 S CONNECT™**

Switch on IMAGE1 S CONNECT™ (see fig. I).

**i NOTE:** The connected LINK modules are turned on automatically.

**7.3.9 Connecting USB storage media**

Connect the USB storage media with the TC 200 IMAGE1 S CONNECT™ to store still images and/or videos (see fig. J).

**i NOTE:** The USB storage medium must have a FAT 32 format.

**7.3.10 Connecting the camera head/ video endoscope**

Insert the camera head connector into the receptacle of the relevant LINK module (see fig. K).

**i NOTE:** Please verify the values set for the primary and secondary sources as well as the layout beforehand (see p. 22).

**i NOTE:** Please note the compatibility list of the respective LINK modules (see chap. 12).

**! WARNING:** Always ensure the connector is completely clean and free of debris. If moisture is present, dry thoroughly with sterile towel prior to insertion.

**7.3.11 Performing the white balance**

**i NOTE:** Perform a white balance before every use. If the same instruments or the same light source were used previously, the white balance may be omitted.

Perform the white balance by pressing the  button on the front of the IMAGE1 S CONNECT™ (TC 200) or by pressing the  in the Live Menu.

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**7.3.12 Test for proper functioning**

1. Inspect the IMAGE1 S module and the camera head and/or the video endoscope for external damage.
2. Ensure that the camera head cable and/or the video endoscope cable is not damaged or kinked.
3. Point the camera head and/or the video endoscope towards an object and check the quality of the display on the monitor.
4. When using 3D video endoscopes, check to make sure the monitor's display is correct. Please ensure that you use the correct polarized eyewear for the monitor.

**i** **NOTE:** *The device and accessories may only be used once the test for proper functioning has been successfully completed.*

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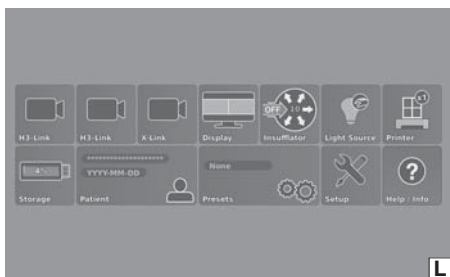
## 8 Operating instructions

### 8.1 Basic information on operation

The use of IMAGE1 S is divided into three menu stages:

- Dashboard (see 8.2)
- Live menu (see 8.3)
- Setup menu (see 8.4)

#### 8.1.1 Keyboard and mouse

The keyboard and mouse can be used to navigate the menus. A keyboard is required to enter patient data.

#### 8.1.2 Camera head buttons

All IMAGE1 FULL HD camera heads feature three buttons for retrieving programmed functions and for menu control:

⑮ scroll up: Scrolls up the menus, options.

⑰ menu/select: Displays the live menu. Once this menu is displayed, this button activates the highlighted option (identical to ENTER on a computer keyboard).

⑯ scroll down: Scrolls down the menus, options.

## 8.2 Dashboard

The Dashboard displays the devices, data storage devices and LINK modules connected to IMAGE1 S CONNECT™ (see fig. L). Unconnected devices are shown in gray.

**NOTE:** The Dashboard is shown when no camera head/video endoscope is connected.

**NOTE:** Depending on the configuration of your system, individual icons may be deactivated or unavailable. The configuration can be modified in the Setup Menu.

**NOTE:** The icons display the current settings.

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
**8.2.1 Description of functions of the symbols used**

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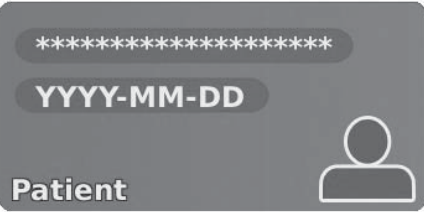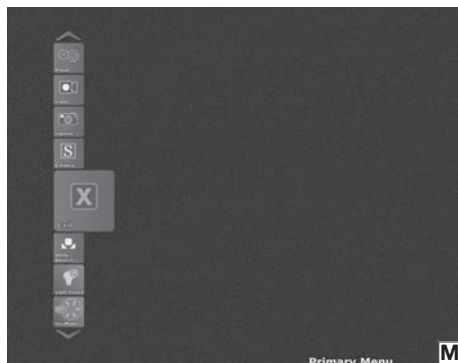
### 8.3 Live Menu

The live menu allows you to choose between different camera functions and change the settings of connected devices (see fig. M). Unconnected devices are shown in gray.

You can access the live menu by pressing the enter button on the keyboard or by pressing the menu button  $\text{M}$  on the camera head/video endoscope or the left side button.

- i** **NOTE:** The icons in the live menu are displayed as soon as a camera head is connected.
- i** **NOTE:** Depending on the configuration of your system, individual icons may be deactivated or unavailable. The configuration can be modified in the Setup Menu.
- i** **NOTE:** The icons display the current settings.

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**8.3.1 Description of functions of the symbols used**

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**i NOTE:** The S technologies are only available with compatible camera heads and video endoscopes (see chapter 12).

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\* SPECTRA A/B is not available in the USA.

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**8.4 Setup Menu**

In the Setup menu, you can change the general settings of the IMAGE1 S CONNECT™ as well as the link-specific settings, and optionally save them as a preset.

You can access the Setup Menu by pressing the setup button on the dashboard as well as by pressing the middle camera head button for 7 seconds in the live menu. Or press Shift-Enter on the keyboard in the setup menu.

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## 9 Cleaning, disinfection, care and sterilization

### 9.1 General safety information

**WARNING:** Risk of infection: Incorrectly reprocessed medical devices expose patients, users and third parties to a risk of infection as well as the risk that the medical device may malfunction. Observe the 'Cleaning, Disinfection, Care, and Sterilization of KARL STORZ Instruments' instructions and the accompanying documentation.

**CAUTION:** When preparing and using the solutions, follow the chemical manufacturer's specifications, paying close attention to proper concentration, exposure time and service life. Prolonged immersion and incorrect concentration may result in damage. Bear in mind the microbiological range of action of the chemicals used.

**CAUTION:** Do not allow water or any liquid to drip inside the unit(s).

**CAUTION:** Danger of damage to medical devices: The use of chemicals which have not been approved by KARL STORZ may cause damage to the medical devices. Only use chemicals approved by KARL STORZ for reprocessing. A complete list can be found on the internet at [www.karlstorz.com](http://www.karlstorz.com).

**CAUTION:** National laws and regulations must be observed.

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**i** **NOTE:** The manual 'Cleaning, Disinfection, Care, and Sterilization of KARL STORZ Instruments' can be downloaded or requested by visiting [www.karlstorz.com](http://www.karlstorz.com).

## 9.2 Wipe-down disinfection of devices

Wipe-clean the exterior surfaces of the medical device with a disposable cloth moistened with disinfectant or a ready-to-use soaked disinfectant cloth. Due to their protein-fixating effect and possible material incompatibility, alcohol-based agents must not be used. The chemical manufacturer's specifications regarding material compatibility must be observed. At the end of the necessary exposure time for the disinfectant as specified by the manufacturer, wipe the surface with a dry low lint cloth.

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### 10.1 Maintenance and safety check

#### 10.1.1 Maintenance

Preventive maintenance is not essential. Regular maintenance can, however, contribute to identifying potential problems before they become serious, thus enhancing the instrument's reliability and extending its useful service life. Maintenance services can be obtained from your local representative or from the manufacturer.

#### 10.1.2 Safety check

Regardless of the accident prevention regulations or testing intervals for medical devices prescribed in different countries, the following checks must be carried out on this device at least once a year by persons who, due to their training, knowledge and practical experience, are able to properly carry out such safety checks and who are not subject to any instructions as far as such checking activities are concerned.

#### Visual inspection:

- Check the device and accessories for any mechanical damage which may impair functionality.
- Check that the labels relevant to safety are legible.

#### Electric measurements:

- Protective grounding measured according to IEC 62353 (current edition), and meets the limits specified in the standard.
- Touch Current, Earth Leakage Current, and Patient Leakage currents, measured according to IEC 62353 (current edition), and meets the limits specified in the standard.

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**Test for proper functioning:**

- A functional check must be carried out in accordance with the instruction manual (see § 7.3.12).
- The safety check must be entered in the unit log-book and the results of the check must be documented.
- If the device does not function reliably and/or safely, it must be repaired.

**10.2 Servicing and repair**

Defective items of equipment must be serviced and repaired exclusively by persons authorized by us; all repair work must employ genuine KARL STORZ parts only.

**10.3 Disposal**

This unit has been marked in accordance with the European Directive on waste electrical and electronic equipment (WEEE).

At the end of its useful operating life, dispose of the unit as electronic scrap.

Please ask either KARL STORZ GmbH & Co. KG, a KARL STORZ subsidiary or your specialist dealer for information on your local collection point.

Within the scope of application of this Directive, KARL STORZ GmbH & Co. KG is responsible for the proper disposal of this device.

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#### 10.4 Repair program

For equipment, individual repair is necessary. To bridge the repair period, you usually receive a unit on loan which you then return as soon as you receive the repaired unit.

In the case of repairs, please contact your local KARL STORZ branch or authorized dealer.

##### **Important information**

In order to prevent infection, it is strictly forbidden to ship contaminated medical devices. All medical devices must be decontaminated on site, to avoid contact and aerogenous infections (among personnel). We reserve the right to return contaminated instruments or devices to the sender.

Repairs, modifications, or expansions which are not performed by KARL STORZ or by experts authorized by KARL STORZ will invalidate all warranty rights. KARL STORZ gives no guarantee on the correct functioning of equipment or instruments which have been repaired by unauthorized third parties.

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### 10.5 Limitation of liability

KARL STORZ shall be liable for failure or deterioration in the safe operation, operational reliability, and performance of this equipment only subject to the following conditions:

- all assembly operations, system expansions, readjustments, modifications or repairs have been performed by a person or persons duly authorized by KARL STORZ;
- the electrical installation of the location in which the device is connected and operated conforms to the applicable statutes and standards, and
- the unit has been used in accordance with its operating instructions at all times.

### 10.6 Warranty

For a period of two years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship. In so doing we cannot accept to bear the cost of transportation or the risk of shipment. The warranty specified in our Standard Terms and Conditions of Business shall additionally apply.

The device is always to be sent to your subsidiary (see chapter 'Subsidiaries') even during the warranty period.

Opening the equipment or performance of any repairs or modifications to the equipment by unauthorized persons shall relieve us of any liability for its performance. Any such opening, repair, or modification performed during the warranty period shall void all warranty.

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### 11.2 Technical documentation

On request, the manufacturer will provide those circuit diagrams, itemized parts listings, descriptions, sets of adjustment instructions and other items of available documentation to suitably qualified user personnel duly authorized by the manufacturer for their use in repairing those components of the unit that have been designated by their respective manufacturers as repairable.

Supply of such technical documentation relating to the unit shall not be construed as constituting manufacturer's authorization of user's personnel, regardless of their level of technical training, to open or repair the unit.

Explicitly exempted herefrom are those maintenance and repair operations described in this manual.

### 11.3 Standard compliance

**According to IEC 60601-1 + A1 + A2, IEC 60601-2-18, UL 60601-1, CSA C22.2 No. 601.1-M90:**

- Type of protection against electric shocks: Protection Class I
- Degree of protection against electric shocks: Applied part of type CF defibrillation-proof 
- Type of protection against moisture: IPX 0

**According to IEC 60601-1-2:**

Please observe the information on electromagnetic compatibility in § 13 (p. 52).

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#### 11.4 Directive compliance for IMAGE1 S (TC 200, TC 300, TC 301 and TC 302)

##### According to Medical Device Directive (MDD):

This medical device belongs to Class I

This medical product bears the CE mark in accordance with the Medical Device Directive (MDD) 93/42/EEC.

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#### 11.5 Software ownership and licensing

With respect to KARL STORZ products containing software components, Customer has a non-exclusive, limited, non-transferable license to use the programmed logic, computer programs and/or software supplied by KARL STORZ, in connection with, and incorporated into, KARL STORZ Products System (Software) internally, only in the form in which delivered to Customer and for the sole purpose of operating in accordance with KARL STORZ's written instructions for such products sold to Customer (and for no other product). The Software, and all modifications, enhancements and upgrades thereto, will, at all times, remain the property of KSEA. Customer may not duplicate, copy, reverse-engineer, de-compile, or disassemble the Software or in any way modify the Software. Customer has no right to, and may not, create derivatives of the software, and Customer may not attempt to copy, create or re-create the source code of the Software. Any and all such modifications or enhancements to the Software by Customer, in contravention of this license, will immediately become the sole property of KSEA.

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Customer hereby acknowledges and agrees (I) that the purchase of KARL STORZ products does not constitute a sale of the Software, (II) that the Software is the property of KARL STORZ, (III) that Customer neither owns nor acquires any claim or intellectual property right in or to the Software, and (IV) KARL STORZ retains all right, title, and interest in and to the Software, at all times, regardless of the form or media in or on which the original or other copies of the Software may exist. In the event of a failure of Customer, or its agents, employees or representatives, to comply with any terms and conditions of the License herein granted, the License will, without any further action by KARL STORZ or any other party, immediately terminate.

**Attribution Notice**

This device contains open source software. You can view a list of the open source software included in the product, and the open source software license terms, on the IMAGE1 S CONNECT™ device itself. From the Dashboard, select “Setup/General”, then “Software Licenses”.

To obtain a copy of the source code being made publicly available by KARL STORZ Imaging, Inc. (“KARL STORZ“) related to open source software used in this KARL STORZ product, you should send your request, along with product part number and serial number, in writing by email to [opensource@karlstorz.com](mailto:opensource@karlstorz.com) or by postal mail to:

KARL STORZ Imaging, Inc.  
Attention: Open Source Software Inquiries  
175 Cremona Drive  
Goleta, CA 93117  
USA

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Name	Article number
IMAGE1 S H3-Z Three-Chip FULL HD Camera Head	TH 100
IMAGE1 S H3-ZI Three-Chip FULL HD Camera Head	TH 101
IMAGE1 S H3-Z FI Three- Chip FULL HD Camera Head	TH 102*
IMAGE1 S H3-P Three-Chip FULL HD Pendelum Camera Head	TH 103
IMAGE1 S H3-ZA Three-Chip FULL HD Camera Head	TH 104
IMAGE1 S H3-M COVIEW® Three-Chip FULL HD Camera Head	TH 106

Naam	Artikelnummer
IMAGE1 S H3-Z 3-chip FULL HD-camerakop	TH 100
IMAGE1 S H3-ZI 3-chip FULL HD-camerakop	TH 101
IMAGE1 S H3-Z FI 3-chip FULL HD-camerakop	TH 102*
IMAGE1 S H3-P 3-chip FULL HD pendelum- camerakop	TH 103
IMAGE1 S H3-ZA 3-chip FULL HD-camerakop	TH 104
IMAGE1 S H3-M COVIEW® 3-chip FULL HD-camerakop	TH 106

**12.2 Kompatible Kameraköpfe  
für IMAGE1 S H3-LINK TC 300  
ohne S-Technologien**
**12.2 Compatible camera heads  
for IMAGE1 S H3-LINK TC 300  
without S technologies**
**12.2 Compatibele camerakoppen  
voor IMAGE1 S H3-LINK  
TC 300 zonder S-technologieën**

Name	Artikelnummer
IMAGE1 H3-Z Drei-Chip FULL HD Kamerakopf	<b>22</b> 2200 55-3
IMAGE1 H3-ZI Drei-Chip FULL HD Inline-Kamerakopf	<b>22</b> 2200 53-3
IMAGE1 H3-Z FI Drei-Chip ICG FULL HD Kamerakopf	<b>22</b> 2200 85-3*
IMAGE1 H3-P Drei-Chip FULL HD Pendelum- Kamerakopf	<b>22</b> 2200 56-3
IMAGE1 H3-ZA Drei-Chip FULL HD Kamerakopf	<b>22</b> 2200 61-3
IMAGE1 H3-FA Drei-Chip FULL HD Kamerakopf	<b>22</b> 2200 60-3

Name	Article number
IMAGE1 H3-Z Three-Chip FULL HD Camera Head	<b>22</b> 2200 55-3
IMAGE1 H3-ZI Three-Chip FULL HD Inline Camera Head	<b>22</b> 2200 53-3
IMAGE1 H3-Z FI Three-Chip ICG FULL HD Camera Head	<b>22</b> 2200 85-3*
IMAGE1 H3-P Three-Chip FULL HD Pendelum Camera Head	<b>22</b> 2200 56-3
IMAGE1 H3-ZA Three-Chip FULL HD Camera Head	<b>22</b> 2200 61-3
IMAGE1 H3-FA Three-Chip FULL HD Camera Head	<b>22</b> 2200 60-3

Naam	Artikelnummer
IMAGE1 H3-Z 3-chip FULL HD-camerakop	<b>22</b> 2200 55-3
IMAGE1 H3-ZI 3-chip FULL HD Inline-camerakop	<b>22</b> 2200 53-3
IMAGE1 H3-Z FI 3-chip ICG FULL HD-camerakop	<b>22</b> 2200 85-3*
IMAGE1 H3-P 3-chip FULL HD pendelum- camerakop	<b>22</b> 2200 56-3
IMAGE1 H3-ZA 3-chip FULL HD-camerakop	<b>22</b> 2200 61-3
IMAGE1 H3-FA 3-chip FULL HD-camerakop	<b>22</b> 2200 60-3

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**12.3 Compatible camera heads  
for IMAGE1 S X-LINK TC 301  
with S technologies**

Name	Article number
IMAGE1 S D1, 50 Hz	TH 115
IMAGE1 S D1, 60 Hz	TH 116

**12.4 Compatible camera heads  
for IMAGE1 S X-LINK TC 301  
without S technologies**

Name	Article number
IMAGE1 HX	TH 110*
IMAGE1 HX-P	TH 111*

**12.5 Compatible video endoscopes  
for IMAGE1 S X-LINK TC 301  
with S technologies**

Name	Article number (PAL/NTSC)
<b>Specialty ENT</b>	
KARL STORZ Video Bronchoscope	11900 BP/BN
Video rhinolaryngoscope	11101 VP/VN
Video rhinolaryngoscope	11101 VPS/VNS
<b>Specialty GASTRO</b>	
Video-Gastroscope 9,3/2,8/1100	13821 PKS/ NKS
Video-Gastroscope Ped./Sten. 5,9/2,0/1100	13820 PKS/ NKS
Video-Gastroscope 12/4,2/1100	13807 PKS/ NKS
Video-Gastroscope 12/2,8 – 3,4/1100	13806 PKS/ NKS
Video-Gastroscope 9,7/2,8/1100	13801 PKS/ NKS
Video-Duodenoscope 12,6/4,2/1260	13885 PKS/ NKS*

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**12.6 Compatible video endoscopes  
and camera heads for  
IMAGE1 S D3-LINK™ TC 302  
with S technology**

Name	Ø [mm]	Article number (PAL/NTSC)
TIPCAM®1 S 3D, 0°	10	26605 AA
TIPCAM®1 S 3D, 30°	10	26605 BA
TIPCAM®1 S 3D, 0°	4	7240 AA3D*
TIPCAM®1 S 3D, 30°	4	7240 BA3D*
VITOM® 3D		TH 200

Naam	Artikelnummer (PAL/NTSC)
<b>Vakgebied PROCTOLOGIE</b>	
Flexibele TROIDL-rectoscoop	13912 PKS/ NKS*
<b>Vakgebied CHIRURGIE</b>	
Videocholedochoscoop	11292 VS/VP
Videocholedochoscoop	11292 VSU/VPU
<b>Vakgebied NEURO</b>	
Video-neuro-endoscoop	11161 V*

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Bestell-Nr. Order no. Bestelnr.	Artikel	Article	Artikel
9826 NB-3D	3D-Monitor, 26"	3D Monitor, 26"	3D-monitor, 26"
9832 NB-3D	3D-Monitor, 32"	3D Monitor, 32"	3D-monitor, 32"
TM 323	3D-Monitor, 32"	3D Monitor, 32"	3D-monitor, 32"

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**13 Information on  
electromagnetic  
compatibility (EMC)**

**WARNING:** Medical electrical devices are subject to special precautions regarding Electromagnetic Compatibility (EMC). Please observe the EMC instructions in this instruction manual during installation and operation.

**WARNING:** IMAGE1 S should not be placed in the immediate vicinity of other electrical equipment. If this cannot be avoided, please ensure that the performance of either product is not affected by unintended electromagnetic coupling.

**WARNING:** IMAGE1 S is only intended for use by physicians and medical assistants. This is a CISPR Class A medical equipment/system. In a domestic environment, IMAGE1 may cause radio interference, in which case it may be necessary to take adequate mitigation measures, such as re-orienting, relocating, or shielding the unit or filtering the connection to the public mains network.

**NOTE:** The EMC tables and other guidelines that are included in this Instruction Manual provide information to the customer or user that is essential in determining the suitability of the Equipment or System for the Electromagnetic Environment of use, and in managing the Electromagnetic Environment of use to permit the Equipment or System to perform its intended use without disturbing other Equipment and Systems or non-medical electrical equipment.

**NOTE:** This device has been tested and complies with the EMC limits values for the Medical Device Directive 93/42/EEC [EN 55011 Class A and IEC 60601-1-2].

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## Information on electromagnetic compatibility (EMC)

These limits are designed to provide reasonable protection against harmful interference in a typical medical installation. The equipment generates uses and can radiate radio frequency energy. If not installed and used in accordance with the instructions, may cause harmful interference to other devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference with other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving device
- Increase the separation between the individual devices
- Connect the equipment into an outlet on a circuit different from that to which the other device(s) is connected
- Consult the manufacturer or field service technician for help.

**WARNING:** The use of portable and mobile RF equipment may have impact on this and other pieces of medical equipment.

**WARNING:** The use of an accessory or cable with the IMAGE1 S other than those specified in this manual may result in increased emissions or decreased immunity of the IMAGE1 S. When using an accessory or cable with IMAGE1 S other than those specified in this manual, it becomes the responsibility of the user of the IMAGE1 S to determine compliance with IEC 60601-1-2 when using this item.

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System Cables and Maximum Lengths used for EMC Compliance				
Type	Shielded	Max Length [m]	Ferrite	Used for
LINK cable	Yes	2	No	Interconnection of inter-module LINK signals
BNC to BNC	Yes	10	No	Connection of 3G-SDI video signals
DVI	Yes	3	No	Connection of DVI signals
4-pin mini-DIN	Yes	1	No	Connection of SCB signals
USB	Yes	2	No	Connection to printer, storage device, keyboard
Power cord	No	3	No	Connection of Device to AC Mains
Ethernet	No	2	No	Connection to LAN (RJ-45) port
IMAGE1 PILOT TC 014	Yes	5	No	Connection to VITOM® 3D

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<b>Table 1</b>		
<b>Guidelines and manufacturer's declaration – electromagnetic emissions</b>		
The IMAGE1 S is intended for use in the electromagnetic environment specified below. The customer or user of the IMAGE1 S should ensure that it is used in such an environment.		
<b>Emissions test</b>	<b>Compliance</b>	<b>Electromagnetic enforcement – guidance</b>
RF emissions acc. to CISPR 11	Group 1	The IMAGE1 S uses RF energy only for its internal function. The customer or user of the IMAGE1 S should make sure that it is used in such an environment. Therefore, its RF emissions are very low and are not likely to cause interference in nearby electronic equipment.
RF Emissions CISPR 11	Class A	The IMAGE1 S is suitable for use in all establishments other than domestic and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

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<b>Table 2</b>			
<b>Guidance and manufacturer's declaration – electromagnetic immunity</b>			
<p>The IMAGE1 S is intended for use in the electromagnetic environment specified below. The user of IMAGE1 S should make sure that it is used in such an environment.</p>			
Immunity tests	EN/IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) acc. to IEC 61000-4-2	± 6 kV Contact discharge ± 8 kV Air discharge	± 6 kV Contact discharge ± 8 kV Air discharge	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Rapid transient electrical interferences/bursts acc. to IEC 61000-4-4	± 2 kV for power lines ± 1 kV for input and output lines	± 2 kV/1 kV ± 1 kV for input and output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surges acc. to IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV Differential mode ± 2 kV Common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and fluctuations of the supply voltage acc. to IEC 61000-4-11	<5 % $U_T^*$ (>95 % dip of the $U_T$ ) for ½ cycle  40 % $U_T$ (60 % dip of the $U_T$ ) for 5 cycles  70 % $U_T$ (30 % dip of the $U_T$ ) for 25 cycles  <5 % $U_T$ (>95 % dip of the $U_T$ ) for 5 seconds	<5 % $U_T^*$ (>95 % dip of the $U_T$ ) for ½ cycle  40 % $U_T$ (60 % dip of the $U_T$ ) for 5 cycles  70 % $U_T$ (30 % dip of the $U_T$ ) for 25 cycles  <5 % $U_T$ (>95 % dip of the $U_T$ ) for 5 seconds	Mains power quality should be that of a typical commercial or hospital environment.  If the user of the IMAGE1 S require continued operation during power mains interruptions, it is recommended that the IMAGE1 S be powered from an uninterruptible power supply or a battery.
Magnetic field at the supply frequency (50/60 Hz) acc. to IEC 61000-4-8	3 A/m	3 A/m	If image distortion occurs, it may be necessary to position the IMAGE1 S further from sources of power frequency magnetic fields or to install magnetic shielding. The power frequency magnetic field should be measured in the intended installation location to assure that it is sufficiently low.
<p> <b>*NOTE:</b> <math>U_T</math> is the mains alternating voltage before applying the test levels.</p>			

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Table 4			
Guidance and manufacturer's declaration – electromagnetic immunity			
The IMAGE1 S is intended for use in the electromagnetic environment specified below. The customer or user of the IMAGE1 S should ensure that it is used in such an environment.			
Immunity tests	EN/IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
			Portable and mobile RF communications equipment should be used no closer to any part of the IMAGE1 S, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance:
Conducted RF disturbances acc. to IEC 61000-4-6	3 V <sub>rms</sub> 150 kHz to < 80 MHz	3 V <sub>rms</sub>	$d = 1.2 \sqrt{P}$
Radiated RF disturbances acc. to IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	$d = 1.2 \sqrt{P}$ 80 MHz to < 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz
			<p>Where P is the nominal power of the transmitter in watts [W] according to the information provided by the manufacturer of the transmitter and d is the recommended separation distance in meters [m].</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,<sup>a</sup> should be less than the compliance level in each frequency range<sup>b</sup>.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<p><b>i</b> <b>NOTE:</b> At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p><b>i</b> <b>NOTE:</b> These guidelines may not apply in all situations. The propagation of electromagnetic quantities is affected by absorptions and reflections of buildings, objects, and persons.</p>			
<p><sup>a</sup> The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz.</p> <p><sup>b</sup> The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2,5 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 has been incorporated into the formulae used in calculating the recommended separation distance for transmitters in these frequency ranges.</p> <p><sup>c</sup> Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the IMAGE1 S is used exceeds the applicable RF compliance level above, the IMAGE1 S should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the IMAGE1 S.</p> <p><sup>d</sup> Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.</p>			

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Table 6 Recommended separation distances between portable and mobile RF communications equipment and the IMAGE1 S			
The IMAGE1 S is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or user of the IMAGE1 S can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the IMAGE1 S as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter [W]	Separation distance d [m] as a function of the transmitter frequency		
	150 kHz to < 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to < 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance D in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
<p> <b>NOTE:</b> At 80 MHz and 800 MHz, the separation distance for higher frequency range applies.</p> <p> <b>NOTE:</b> These guidelines may not apply in all situations. The propagation of electromagnetic quantities is affected by absorptions and reflections of structures, objects and people.</p>			

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KARL STORZ Endoscopy Canada Ltd.  
7171 Millcreek Drive, Mississauga  
Ontario, L5N 3R3, Canada  
Phone: +1 905 816-4500, Fax: +1 905 858-4599  
Toll free phone: 1-800-268-4880 (Canada only)  
Toll free fax: 1-800-482-4198 (Canada only)  
E-Mail: info.canada@karlstorz.com

KARL STORZ Endoscopy-America, Inc.  
2151 East Grand Avenue  
El Segundo, CA 90245-5017, USA  
Phone: +1 424 218-8100, Fax: +1 424 218-8526  
Toll free phone: 800 421-0837 (USA only)  
Toll free fax: 800 321-1304 (USA only)  
E-Mail: info@kseas.com

KARL STORZ Veterinary Endoscopy-America, Inc.  
1 South Los Cameros Road  
Goleta, CA 93117, USA  
Phone: +1 805 968-7776, Fax: +1 805 685-2588  
E-Mail: info@karlstorzvet.com

KARL STORZ Endoscopia Latino-America, Inc.  
815 N. W. 57th Avenue, Suite 480  
Miami, FL 33126-2042, USA  
Phone: +1 305 262-8980, Fax: +1 305 262-8986  
E-Mail: info@ksela.com

KARL STORZ Endoscopia  
Miramar Trade Center,  
Edificio Jerusalem, Oficina 108,  
La Habana, Cuba  
Phone: +53 72041097, Fax: +53 72041098  
E-Mail: karlstorz@enet.cu

KARL STORZ Endoscopia México S.A de C.V.  
Lago Constanza No 326,  
Col. Granada Del. Miguel Hidalgo  
C.P. 11520 México D.F.  
Phone: +52 (55) 1101 1520  
E-Mail: mx-info@karlstorz.com

KARL STORZ Marketing América Do Sul Ltda.  
Rua Joaquim Floriano, nº. 413, 20º andar – Itaim Bibi,  
CEP-04534-011 São Paulo, Brasil  
Phone: +55 11 3526-4600, Fax: +55 11 3526-4680  
E-Mail: info@karlstorz.com.br

KARL STORZ Endoscopia Argentina S.A.  
Zufriategui 627 6° Piso, B1638 CAA - Vicente Lopez  
Provincia de Buenos Aires, Argentina  
Phone: +54 11 4718 0919, Fax: +54 11 4718 2773  
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskopi Norge AS  
Østensjøveien 15B  
0661 Oslo, Norway  
Phone: +47 6380 5600, Fax: +47 6380 5601  
post@karlstorz.no

KARL STORZ Endoskop Sverige AB  
Storsåtragränd 14  
12739 Skärholmen, Sweden  
14108 Kungens Kurva, Sweden  
Phone: +46 8 505 648 00  
E-Mail: kundservice@karlstorz.se

KARL STORZ Endoscopy Suomi OY  
Valimotie 1 B, 4. kerros  
00380 Helsinki, Finland  
Phone: +358 (0)96824774, Fax: +358 (0)968247755  
E-Mail: asiakaspalvelu@karlstorz.fi

KARL STORZ GmbH & Co. KG  
Representation Office  
Kestucio g. 59  
08124 Vilnius, Lithuania  
Phone: +370 68567000  
E-Mail: ausra.kukeniene@karlstorz.com

KARL STORZ Endoskopi Danmark A/S  
Skovlyoftens 33  
2840 Holte, Denmark  
Phone: +45 45162600, Fax: +45 45162609  
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.  
415 Perth Avenue, Slough  
Berkshire, SL1 4TQ, United Kingdom  
Phone: +44 1753 503500, Fax: +44 1753 578124  
E-Mail: info-uk@karlstorz.com

KARL STORZ Endoscopie Nederland B. V.  
Displayweg 2  
3821 BT Amersfoort, Netherlands  
Phone: +31 (0)33 4545890  
E-Mail: info-nl@karlstorz.com

KARL STORZ Endoscopy Belgium N. V.  
Phone: +31 (0)33 4545890  
E-Mail: info-be@karlstorz.com

KARL STORZ Endoscopie France S. A. S.  
12, rue Georges Guynemer, Quartier de l'Europe  
78280 Guyancourt, France  
Phone: +33 1 30484200, Fax: +33 1 30484201  
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH  
Landstraßer Hauptstr. 148/1/G1  
1030 Wien, Austria  
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479  
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Ibérica S. A.  
Parque Empresarial San Fernando  
Edificio Munich – Planta Baja  
28830 Madrid, Spain  
Phone: +34 91 6771051, Fax: +34 91 6772981  
E-Mail: info-es@karlstorz.com

KARL STORZ Endoscopia Italia S. r. l.  
Via dell'Artigianato, 3  
37135 Verona, Italy  
Phone: +39 045 8222000, Fax: +39 045 8222001  
E-Mail: info-ita@karlstorz.com

KARL STORZ Adria Eos d.o.o.  
Capraška 6  
10000 Zagreb, Croatia  
Phone: +385 1 6406 070, Fax: +385 1 6406 077  
E-Mail: info@karlstorz.hr

KARL STORZ Endoskopija d.o.o.  
Verovškova c. 60A  
1000 Ljubljana, Slovenia  
Phone: +386 1 620 5880, Fax: + 386 1 620 5882  
E-Mail: pisarna@karlstorz.si

KARL STORZ Endoskop Polska Marketing Sp. z o.o.  
ul. Bojkowska 47  
44-100 Gliwice, Poland  
Phone: +48 32 706 13 00, Fax: +48 32 706 13 07  
E-Mail: info-pl@karlstorz.com

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KARL STORZ ENDOSKOPE  
East Mediterranean and Gulf S.A.L.  
Block M, 3rd Floor  
Beirut Souks, Weygand Street  
2012 3301 Beirut, Lebanon  
Phone: +961 1 999390, Fax +961 1 999391  
E-Mail: info@karlstorz-emg.com

KARL STORZ Endoscopy (South Africa) (Pty) Ltd.  
P.O. 6061, Roggebaai 8012  
Cape Town, South Africa  
Phone: +27 21 417 2600, Fax: +27 21 421 5103  
E-Mail: info@karlstorz.co.za

TOO KARL STORZ Endoscopy Kasachstan  
6. Saryarka str. BC „Arman“, office 314  
010000 Astana, Republic of Kazakhstan  
Phone: +7 7172 552-549, 552-788, Fax: -444  
E-Mail: info@karlstorz.kz

KARL STORZ ENDOSKOPE  
East Mediterranean & Gulf (branch)  
Building West Side 7A – Unit 7WA – 3008  
Dubai Airport Free Zone, P.O. Box 54983  
Dubai - United Arab Emirates  
Phone: +971 (0)4 2958887, Fax: +971 (0)4 3205282  
Service Hotline: +971 (0) 3415882  
E-Mail: info-gne@karlstorz-emg.com

KARL STORZ Endoscopy India Private Limited  
11th Floor, Dr. Gopal Das Bhawan  
28, Barakhamba Road  
New Delhi 110001, India  
Phone: +91 11 4374 3000, Fax: +91 11 4374 3010  
E-Mail: corporate@karlstorz.in

KARL STORZ GmbH & Co. KG  
Resident Representative Office  
14th Floor, Kumho Asiana Plaza  
39 Le Duan, District 1  
Ho Chi Minh City, Vietnam  
Phone: +84 8 3823 8000, Fax: +84 8 3823 8039  
E-Mail: infovietnam@karlstorz.com

KARL STORZ Endoscopy Australia Pty. Ltd.  
15 Orion Road Lane Cove NSW 2066  
P O Box 50 Lane Cove NSW 1595, Australia  
Phone: +61 (0)2 9490 6700  
Toll free: 1800 996 562 (Australia only)  
Fax: +61 (0)2 9420 0695  
E-Mail: info@karlstorz.au

KARL STORZ Endoscopy China Ltd., Hong Kong  
Room 2503-05, 25F AXA Tower, Landmark East,  
No. 100 How Ming Street, Kwun Tong, Kowloon,  
Hong Kong, People's Republic of China  
Phone: +852 28 65 2411, Fax: +852 28 65 4114  
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy (Shanghai) Ltd.  
Beijing Branch,  
Room 1805-1807, Building B, 18F Beijing IFC  
No. 8, Jianguomenwai Street, Chaoyang District,  
Beijing, 100022, People's Republic of China  
Phone: +86 10 5638188, Fax: +86 10 5638199  
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.  
Shanghai Branch  
36F, Tower 1 Grand Gateway  
No. 1 Hongguiao Road, Xuhui District,  
Shanghai, 200030, People's Republic of China  
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199  
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.  
Chengdu Branch  
Room 805, 8F Jin Jiang Intenational Building  
No. 1 West Linjiang Road, Wuhou District,  
Chengdu, 6100414, People's Republic of China  
Phone: +86 28 86587977, Fax: +86 28 86587975  
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.  
Shenyang Branch,  
Room 2001-2005, 20F N-MEDIA International Center,  
No. 167 Youth Avenue, Shenhe District,  
Shenyang, 110014, People's Republic of China  
Phone: +86 24 23181118, Fax: +86 24 23181119  
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy (Shanghai) Ltd.  
Guangzhou Branch  
Room 02B, 35F Teem Tower,  
No. 208 Tianhe Road, Tianhe District,  
Guangzhou, 510620, People's Republic of China  
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286  
E-Mail: info@karlstorz.com.cn

KARL STORZ Endoscopy Asia Marketing Pte Ltd.  
No. 8 Commonwealth Lane #03-02  
Singapore 149555, Singapore  
Phone: +65 69229150, Fax: +65 69229155  
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd  
No. 8 Commonwealth Lane #03-02  
Singapore 149555, Singapore  
Phone: +65 69229150, Fax: +65 69229155  
E-Mail: infoasia@karlstorz.com.sg

KARL STORZ Endoscopy Korea Co. Ltd.  
I-Park Bundang1, 102dong, 8F  
239 Jeongjail-Ro, Bundang-Gu  
Seongnam, 463-859, Korea  
Phone: +82-70-4350-7474, Fax: +82-70-8277-3299  
E-Mail: stephen.Lee@karlstorz.com

KARL STORZ Endoscopy Taiwan Ltd.  
12F, No. 192, Sec. 2, Chung Hsin Rd.,  
Sindian District, New Taipei City, Taiwan  
Phone: +886 933 014 160, Fax: +886 2 8672 6399  
E-Mail: mingwang@karlstorz.tw

KARL STORZ Endoscopy Japan K. K.  
Stage Bldg. 8F, 2-7-2 Fujimi  
Chiyoda-ku, Tokyo 102-0071, Japan  
Phone: +81 3 6380-8622, Fax: +81 3 6380-8633  
E-Mail: info@karlstorz.jp

www.karlstorz.com

**STORZ**  
KARL STORZ—ENDOSKOPE

 **KARL STORZ GmbH & Co. KG**

Mittelstraße 8  
78532 Tuttlingen  
Postfach 230  
78503 Tuttlingen  
Germany

Telefon: +49 (0)7461 708-0  
Telefax: +49 (0)7461 708-105  
E-Mail: [info@karlstorz.com](mailto:info@karlstorz.com)  
Web: [www.karlstorz.com](http://www.karlstorz.com)