

...Recovery and Patient transport tables - not height adjustable

Please choose - without surcharge – between 2 different variants of leatherettes:
- **Stamskin®Top (34 colours)** or - **Skai® Pandora Plus (20 colours)**

As our standard: Extra upholstery for Recovery and Patient transport tables

The upholstery of this tables is a special combination of a solid foam base coated with a high resilience (cold) foam. Therefore it is more flexible than upholsteries of the standard examination tables. (Thickness of this upholstery: 8,5 cm). To get more comfort for the patients (or lying times more than 2 hours) we recommend a thicker and more resilience upholstery in cold-foam with restoring force, thickness 10 cm.



Picture with extras: Storage basket (Mod. 1850)

Model	Width Lying surface	Width overall: with side guards	Frame colour
18-01	650 mm	750 mm	Chromed
18-00	650 mm	750 mm	White powder painted
18-04	650 mm	750 mm	Silver powder painted
19-01	800 mm	900 mm	Chromed
19-00	800 mm	900 mm	White powder painted
19-04	800 mm	900 mm	Silver powder painted

Model 18...and Model 19....

Underpart of square tube 30x30x2 mm.

Basic equipment:

Lying surface 2-sectioned:

- Head part upwards adjustable (+30°) by 2 metal ratchets, safety and for longevity. Length of head part: 550 mm
- Foot part fixed. Length: 1.400 mm

- With transverse and longitudinal bracing

Longitudinal bracing as additional stabilizer and also as support for holder for storage basket (= option: Mod.-No. 1850)

- Movable on 4 castors, 2 with total lock in diagonal placement.

(Different possibilities for mobility see list of extras, page 34)

- 1 push handle, chromed, standard placed at foot end of the table (possible option: placed at head part or 2 push handles)

- **2 side guards stickable, chromed.**

Length side guard: 700 mm, Height over upholstery: 200 mm.

(Different possibilities for other options of side guards see list of extras, page 34)

Measures: Length: 1.950 mm,

Length overall (inclusive of 1 push handle): 2.080 mm

Height: 680 mm (other heights possible on request)

Load capacity: 250 kg (max. weight of patient)



Side guard Model 1005

to fit on side rail

including attachment clamp.

Side guard could be lowered or moved along side rail.

Measures: L: 800 mm, H: 220 mm

(Side rail not included)

As option: Deliverable in customize measures (For example: Table is accommodated to fit for elevator)



Both pictures with extras: Comfort-casters 125 mm diam., with central lock, castors antistatic (Mod. 1186)

Model 80...and 81....

Underpart made of square tube.

Basic equipment:

Lying surface 2-sectioned:

- Head part stepless upwards adjustable by gas spring till +75° as a back for seat positioning. Length: 750 mm
- Foot part fixed. Length: 1.200 mm

- With transverse and longitudinal bracing

- 1 push handle, chromed, standard placed at foot end of the table (possible option: placed at head part or 2 push handles)

- **2 side guards foldable, chromed.**

Length side guard: 1.060 mm, Height over upholstery: 220 mm

(The downwards folded side guards lay down besides the upholstery, below upholstery surface level. The side guards are fixed in folded up- or downwards position by a click-stop device)

- Movable on comfort-casters 125 mm diam., with central lock, 1 castor with direction lock. For an easier operation of the central lock system there are 4 levers at the underpart, one at each foot.


- 4 Bumpers at the top frame

Measures: Length: 1.950 mm, Length overall: 2.080 mm

Height: 680 mm (other heights possible on request)

Load capacity: 250 kg (max. weight of patient)

Model	Width Lying surface	Width overall: with side guards	Frame colour
80-00	650 mm	800 mm	White powder painted
80-04	650 mm	800 mm	Silver powder painted
81-00	800 mm	950 mm	White powder painted
81-04	800 mm	950 mm	Silver powder painted

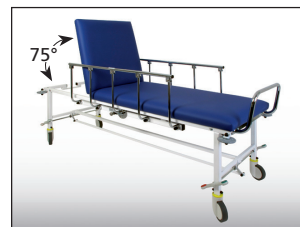

Model 80-00 with extras

For customization of our Recovery and Patient transport tables with special equipment please have a look at page 34 of this catalogue.
For example: - Other side guard systems, - Different possibilities of castors (sizes or locking mechanism),
- Side rails, - Infusion pole, - Upholstery in soft and elastic configuration, - Bumpers, - Paper roll holder.

...Recovery and Patient transport tables - not height adjustable

Please choose - without surcharge – between 2 different variants of leatherettes:
- **Stamskin®Top (34 colours)** or - **Skai® Pandora Plus (20 colours)**

As our standard: Extra upholstery for Recovery and Patient transport tables

The upholstery of this tables is a special combination of a solid foam base coated with a high resilience (cold) foam. Therefore it is more flexible than upholsteries of the standard examination tables. (Thickness of this upholstery: 8,5 cm). To get more comfort for the patients (or lying times more than 2 hours) we recommend a thicker and more resilience upholstery in in cold-foam with restoring force, thickness 10 cm.



Model 118-00

Picture with extras: - 2 side guards up/downwards tiltable, L: 1 m (2x Mod. 1044), - Paper roll holder (Mod. 1210), Wall protection bumpers (4x Mod. 1170), - 2 side rails, L: 350 mm (2x Mod. 1410), 2 special push handles higher configuration

Model 118-01: Frame chromed
Model 118-00: Frame white powder coated
Model 118-04: Frame silver powder coated
Underpart of square tube 30x30x2 mm.

Basic equipment:

Lying surface 2-sectioned:

- Head part upwards adjustable (+30°) by 2 metal ratchets, safety and for longevity. Length of head part: 550 mm
- Foot part fixed. Length: 1.400 mm
- With transverse and longitudinal bracing
- With storage basket
- Movable on 4 castors, 2 with total lock in diagonal placement. (Other possibilities for mobility see list of extras, page 34)
- 2 push handles, chromed (picture shows special push handles in higher configuration, no standard, possible on request)
- **2 side guards stickable**, chromed (L: 700 mm, H: over upholstery: 200 mm) (Different possibilities for other side guards see list of extras, page 34. Picture shows side guards up/downwards tiltable, length 1 m)

Measures lying surface: L: 1.950 mm, W: 650 mm, H: 650 mm

Measures overall: L: 2.200 mm, W: 750 mm

Load capacity: 250 kg (max. weight of patient)



Model 180-00

Picture with extras: - Wall protection bumpers (4x Mod. 1170)

Model 180-00: Frame white powder coated
Model 180-04: Frame silver powder coated
Underpart of square tube 30x30x2 mm.

Basic equipment:

Lying surface 3-sectioned:

- Head part stepless upwards adjustable till +75° by gas spring. Length head part: 750 mm.
- Middle part (Length: 300 mm) not adjustable
- Foot part stepless upwards adjustable by gas spring till +40°. Length foot part: 900 mm
- With 4 corner bracings for additional stabilization (Therefore free space under table top)
- Movable on 4 castors, each with total lock. (Other possibilities for mobility see list of extras, page 34. Please note, castors with central lock are not possible)
- 1 push handle, chromed, standard placed at foot end
- **2 side guards up/downwards tiltable**

Length: 700 mm, Height over upholstery: 200 mm

Measures lying surface: L: 1.950 mm, W: 650 mm, H: 650 mm

Measures overall: L: 2.080 mm, W: 750 mm

Load capacity: 250 kg (max. weight of patient)



Model 188-01

Picture with extras: - 1 push handle in upright design (Mod. 1815)

Model 188-01: Frame chromed
Model 188-00: Frame white powder coated
Model 188-04: Frame silver powder coated
Underpart of square tube 30x30x2 mm.

Basic equipment:

Lying surface 2-sectioned:

- Head part (L: 550 mm) upwards adjustable +30° by 2 metal ratchets
- Foot part fixed. Length: 1.400 mm
- With 4 corner bracings for additional stabilization (Therefore free space under table top)
- Movable on 4 castors, each with total lock. Please note, castors with central lock are not possible) (Other possibilities for mobility see list of extras, page 34.
- 1 push handle, chromed, standard placed at foot end (Note: Picture shows push handle in upright design. This is interesting if the overall measures of the table must be accommodated, for example for elevator.)

- **2 side guards foldable**, chromed

Length side guard: 1.060 mm, Height over upholstery: 220 mm

Measures lying surface: L: 1.950 mm, W: 650 mm,


H: 650 mm (other heights possible)

Measures overall: L: 2.080 mm, W: 770 mm

Load capacity: 250 kg (max. weight of patient)

All tables are produced in our company in Bielefeld, Germany. On request the tables may be customized: Special measures overall to accommodate the table to the local place, special heights, push handles in special ergonomic design, (to have ergonomic push handle height while moving the table), additional upholstery segments and functions. **Please contact us.**

...Recovery and Patient transport tables - height adjustable

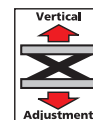

Please choose - without surcharge – between 2 different variants of leatherettes:
- **Stamskin®Top (34 colours)** or - **Skai® Pandora Plus (20 colours)**

As our standard: Extra upholstery for Recovery and Patient transport tables

The upholstery of this tables is a special combination of a solid foam base coated with a high resilience (cold) foam. Therefore it is more flexible than upholsteries of the standard examination tables. (Thickness of this upholstery: 8,5 cm). To get more comfort for the patients (or lying times more than 2 hours) we recommend a thicker and more resilience upholstery in cold-foam with restoring force, thickness 10 cm.



Model 2120-00XL/H

Picture with extras: - Comfort-Castors 125 mm diam. with central locking system (Mod. 1185)
- Storage basket (Mod. 1850)

Model 2120-...XL and Model 2121-...XL

Alternative 1: Manual foot operated height adjustment by hydraulic pump.

Pedals for foot-controlled height positioning are bothsides the table.

Alternative 2: Height adjustment by electromotive drive.

The electric motor is controlled by 1 foot switch

(possible option: Hand switch)

Time of electric level positioning from min. to max.: 22 sec.

Basic equipment:

Lying surface 2-sectioned:

- Head part upwards adjustable +30° by 2 metal rachets, safety and for longevity. Length of head part: 550 mm

- Foot part fixed. Length: 1.400 mm

- 2 push handles, chromed

- **2 side guards** up/downwards **tiltable**

Length side guard: 700 mm, Height over upholstery: 200 mm

(Other side guards, for example side guards foldable, possible)

- Movable on 4 castors, each with total lock

(Other possibilities for mobility see for example picture or page 34)

- 4 Bumpers 150 mm diam. at the edges of top frame

Length lying surface: 1.950 mm, **Length overall:** 2.200 mm

Load capacity: 250 kg (max. weight of patient)



Model 2118-00XL/H

Model 2118-...XL and Model 2119-...XL

like Mod. 2120-...XL and Model 2121-...XL, but difference:

With: - **2 side guards stickable**, chromed.

Length: 700 mm, Height over upholstery: 200 mm.

If the side guards are not used, they could be stucked in converse (showing downwards) into the holders.

Length lying surface: 1.950 mm, **Length overall:** 2.200 mm

Load capacity: 250 kg (max. weight of patient)

For customization of our Recovery and Patient transport tables with special equipment please have a look at page 34 of this catalogue. For example:

- Other side guard systems, - Different possibilities of castors (different sizes and locking mechanism),
- Side rails, - Infusion pole, - Upholstery in soft and elastic configuration, - Bumpers, - Storage basket,
- Paper roll holder, Head part stepless adjustable by gas spring
- Special measures overall to accommodate the table to the local place


Please do not hesitate to contact us.

Models with side guards stickable	Models with side guards up/downwards tiltable	Width Lying surface	Width overall: With side guards	Vertical height adjustment by:	Frame colour
2118-00XL/H	2120-00XL/H	650 mm	750 mm	Hydraulic pump: 520 - 940 mm	White powder painted
2118-04XL/H	2120-04XL/H	650 mm	750 mm	Hydraulic pump: 520 - 940 mm	Silver powder painted
2119-00XL/H	2121-00XL/H	800 mm	900 mm	Hydraulic pump: 520 - 940 mm	White powder painted
2119-04XL/H	2121-04XL/H	800 mm	900 mm	Hydraulic pump: 520 - 940 mm	Silver powder painted
2118-00XL	2120-00XL	650 mm	750 mm	Electric motor: 500 - 940 mm	White powder painted
2118-04XL	2120-04XL	650 mm	750 mm	Electric motor: 500 - 940 mm	Silver powder painted
2119-00XL	2121-00XL	800 mm	900 mm	Electric motor: 500 - 940 mm	White powder painted
2119-04XL	2121-04XL	800 mm	900 mm	Electric motor: 500 - 940 mm	Silver powder painted

Tables with electric motor powered adjustment must have a safety-lock-box to protect against unauthorized usage (regulations of BfArM based on EU-ordinance acc. EN 60601). Please choose between our 2 possibilities:

- Mod. 1993: Safety box, fixed on the underpart, with 1 magnetic key (see picture on page 34)
- Mod. 1995: Emergency stop button, fixed on top frame, with 2 keys (see picture on page 34)

...Recovery and Patient transport tables - height adjustable

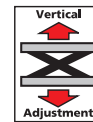

Please choose - without surcharge – between 2 different variants of leatherettes:
- **Stamskin®Top (34 colours)** or - **Skai® Pandoria Plus (20 colours)**

As our standard: Extra upholstery for Recovery and Patient transport tables

The upholstery of this tables is a special combination of a solid foam base coated with a high resilience (cold) foam. Therefore it is more flexible than upholsteries of the standard examination tables. (Thickness of this upholstery: 8,5 cm). To get more comfort for the patients (or lying times more than 2 hours) we recommend a thicker and more resilience upholstery in cold-foam with restoring force, thickness 10 cm.



Model 2210-00XL/H

Picture with extras:

- One-part foot-operating lever for central locking (Mod. 1199)
- Storage basket (Mod. 1850)
- 2 side rails, each 350 mm length (2x Mod. 1410)

Model 2210-...XL and Model 2215-...XL

Manual foot operated height adjustment by hydraulic pump from 630 mm to 1.050 mm.

Pedals for foot-controlled height positioning are bothsides the table.

Basic equipment:

Lying surface 2-sectioned:

- Head part stepless upwards adjustable by gas spring till +75°, as a back for seat positioning.

Length head part: 750 mm

- Foot part fixed.

Length foot part: 1.200 mm

- 2 push handles, chromed

- **2 side guards foldable**, chromed.

Length side guard: 1.060 mm, Height over upholstery: 220 mm

(The downwards folded side guards lay down besides the upholstery, below upholstery surface level. The side guards are fixed in folded up- or downwards position by a click-stop device)

- Movable on comfort-casters with double wheels 125 mm diam., with central lock, 1 castor with direction lock.

For an easier operation of the central lock system there are 4 levers at the underpart, one at each foot.

- 4 Bumpers 200 mm diam. at the edges of top frame

Length lying surface: 1.950 mm

Length overall: 2.200 mm

Load capacity: 250 kg (max. weight of patient)



Model 2250-00XL/H

Picture with extras: - Paper roll holder (Mod. 1267)

Model 2250-...XL/H and Model 2255-...XL/H

like Mod. 2210-...XL/H and Mod. 2215-...XL/H, but difference:

Lying surface 3-sectioned:

- Head part stepless upwards adjustable by gas spring till +75°, as a back for seat positioning.

Length head part: 750 mm

- Middle part fixed.

Length middle part: 300 mm

- Foot part stepless upwards adjustable by gas spring till +35°.

Length foot part: 900 mm

Height adjustment: 630 mm to 1.050 mm

Length lying surface: 1.950 mm

Length overall: 2.200 mm

Load capacity: 250 kg (max. weight of patient)

All tables are produced in our company in Bielefeld, Germany. On request the tables may be deliverable in special measures (Length, Width, Height) or with additional functions. Please also have a look at our special equipment on page 34 in this catalogue to provide the tables for your special applications.

Please do not hesitate to contact us.


Models with lying surface 2-sectioned	Models with lying surface 3-sectioned	Width lying surface	Width overall: with side guards	Height adjustment vertical by:	Frame colour
2210-00XL/H	2250-00XL/H	650 mm	800 mm	Hydraulic pump: 630 - 1.050 mm	White powder painted
2210-04XL/H	2250-04XL/H	650 mm	800 mm	Hydraulic pump: 630 - 1.050 mm	Silver powder painted
2215-00XL/H	2255-00XL/H	800 mm	950 mm	Hydraulic pump: 630 - 1.050 mm	White powder painted
2215-04XL/H	2255-04XL/H	800 mm	950 mm	Hydraulic pump: 630 - 1.050 mm	Silver powder painted

Height adjustment by electric motor drive possible on request

...Recovery and Patient transport tables

- height adjustable

- Trendelenburg adjustment of the top

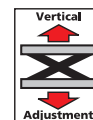

As our standard: Extra upholstery for Recovery and Patient transport tables

The upholstery of this tables is a special combination of a solid foam base coated with a high resilience (cold) foam. Therefore it is more flexible than upholsteries of the standard examination tables. (Thickness of this upholstery: 8,5 cm). To get more comfort for the patients (or lying times more than 2 hours) we recommend a thicker and more resilience upholstery in cold-foam with restoring force, thickness 10 cm.



Model 2420-00XL/H

Model 2420-...XL and Model 2421-...XL

Manual foot operated height adjustment by hydraulic pump from 570 mm to 990 mm.
Pedals for foot-controlled height positioning are bothsides the table.

Basic equipment:

Lying surface 2-sectioned:

- Head part (L: 550 mm) stepless upwards adjustable by gas spring till +60°.
- Foot part fixed. Length: 1.400 mm
- Lying surface adjustable by gas spring to inclined position (Trendelenburg) till 12°
- 2 push handles, chromed

- 2 side guards up/downwards **tiltable**, chromed.

Length side guard: 700 mm, Height over upholstery: 200 mm (Other side guards, for example side guards foldable, are possible)

- Movable on 4 castors 100 mm diam., each with total lock (Other possibilities for mobility, please see page 34)

- 4 Bumpers 150 mm diam., at the edges of top frame

Length lying surface: 1.950 mm. **Length overall:** 2.200 mm

Load capacity: 250 kg (max. weight of patient)

Model 2418-...XL and Model 2419-...XL (without picture)

like Mod. 2420...XL . and 2421...XL, but difference:

With - 2 side guards stickable, chromed.



Both picture with extras:
Castors 100 mm diam. with central lock system Mod.1180N

Model 2420-00XL/H

Models with side guards stickable	Models with side guards up/downwards tiltable	Width Lying surface	Width overall: with side guards	Frame colour
2418-00XL/H	2420-00XL/H	650 mm	750 mm	White powder painted
2418-04XL/H	2420-04XL/H	650 mm	750 mm	Silver powder painted
2419-00XL/H	2421-00XL/H	800 mm	900 mm	White powder painted
2419-04XL/H	2421-04XL/H	800 mm	900 mm	Silver powder painted

Height adjustment by electric motor drive possible on request



Model 2448-00XL/H

Model 2448-...XL and Model 2449-...XL

Manual foot operated height adjustment by hydraulic pump from 570 mm to 990 mm.

Pedals for foot-controlled height positioning are bothsides the table.

Basic equipment:

Lying surface 2-sectioned:

- Head part stepless upwards adjustable by gas spring till +75°, as a back for seat positioning. Length: 750 mm
- Foot part fixed. Length: 1.200 mm

- Lying surface adjustable by gas spring to inclined position (Trendelenburg) till 12°

- 2 push handles, chromed

- 2 side guards stickable, chromed

Length side guard: 700 mm, Height over upholstery: 200 mm.

(Other side guards, for example side guards foldable, are possible)

- Movable on 4 castors 100 mm diam., each with total lock

(Other possibilities for mobility, please see page 34)

- 4 Bumpers 150 mm diam., at the edges of top frame

Length lying surface: 1.950 mm. **Length overall:** 2.200 mm

Load capacity: 250 kg (max. weight of patient)

Model 2450-...XL and Model 2451-...XL (without picture)

like Mod. 2448...XL and Mod. 2449...XL, but difference:

With: - 2 side guards up/downwards tiltable, chromed



Both picture with extras:
Castors 100 mm diam. with central lock system Mod.1180N

Model 2448-00XL/H

Models with side guards stickable	Models with side guards up/downwards tiltable	Width Lying surface	Width overall: with side guards	Frame colour
2448-00XL/H	2450-00XL/H	650 mm	750 mm	White powder painted
2448-04XL/H	2450-04XL/H	650 mm	750 mm	Silver powder painted
2449-00XL/H	2451-00XL/H	800 mm	900 mm	White powder painted
2449-04XL/H	2451-04XL/H	800 mm	900 mm	Silver powder painted

Height adjustment by electric motor drive possible on request

To optimize the handling of Models 2418XL to 2451XL we recommend:
Movability on Comfort-Castors diam. 125 or diam. 150 mm, with central lock.
Please notice: Because of technical reasons the measures of height adjustment changes.


...Recovery and Patient transport tables

- height adjustable

- Trendelenburg adjustment of the top

Please choose - without surcharge – between 2 different variants of leatherettes:

- Stamskin®Top (34 colours) or - Skai® Pandoria Plus (20 colours)

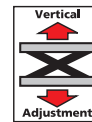

Model 2410-00XL/H

Both picture with extras:
- Storage basket (Mod. 1850)



Model 2410-00XL/H

Model 2410-00XL/H: Frame white powder coated

Model 2410-04XL/H: Frame silver powder coated

Manual foot operated height adjustment by hydraulic pump from **690 mm to 1.110 mm**.

Pedals for foot-controlled height positioning are both sides the table.

Lying-surface with one-piece cover:

- Head part stepless upwards adjustable by gas spring till +70°, as a back for seat positioning.

Length: 720 mm

- Foot part not adjustable.

Length: 1.200 mm

- Lying surface adjustable by gas spring to inclined position (Trendelenburg) till 12°

- With chromed push handles at head and foot end:

Push handle at head end 2 parts, each assembled outside for better access to the patient

Push handle at foot end 1-part over complete width of the table

- 2 side guards foldable, chromed.

Length side guard: 1.060 mm, Height over upholstery: 200 mm

(The downwards folded side guards lay down besides the upholstery, below upholstery surface level. The side guards are fixed in folded up- or downwards position by a click-stop device)

- Movable on comfort-castors 150 mm diam., with central lock, 1 castor with direction lock.

(Castors one-wheeled, with precision-ball-bearings)

- All, 4 castors, are antistatic

- Central and directional locking of castors by 2 change levers at subframe (head end) and one-part operating foot-lever, chromed (foot end)

- 4 bumpers 200 mm diam. at the edges of top frame

- **With detachable mattress, 10 cm thick of high quality resilience (cold) foam**

Mattress with Velcro fastening at the underside.

- Support plate for mattress made of Trespa, white

Measures lying surface: Length: 1.930 mm, Width: 650 mm

Height adjustable from 680 mm to 1.010 mm

Measures overall: Length: 2.030 mm, Width: 800 mm

Load capacity: 250 kg (max. weight of patient)

Some possibilities of accessories as special equipment:

- Side rails, each 350 mm long, bothsides the head part (2 x Mod. 1410)
- Fixing clamp for side rail (Mod. 1430)
- Infusion pole (Mod. 1450)
- Storage basket metal-made, powder painted in colour of underpart frame (Mod. 1850)

All Dewert-Tables are manufactured in Bielefeld, Germany. If possible, the tables may be delivered in special measures or with additional functions for your special applications. **Please do not hesitate to contact us.**